
LED facial mask

en Instructions for use

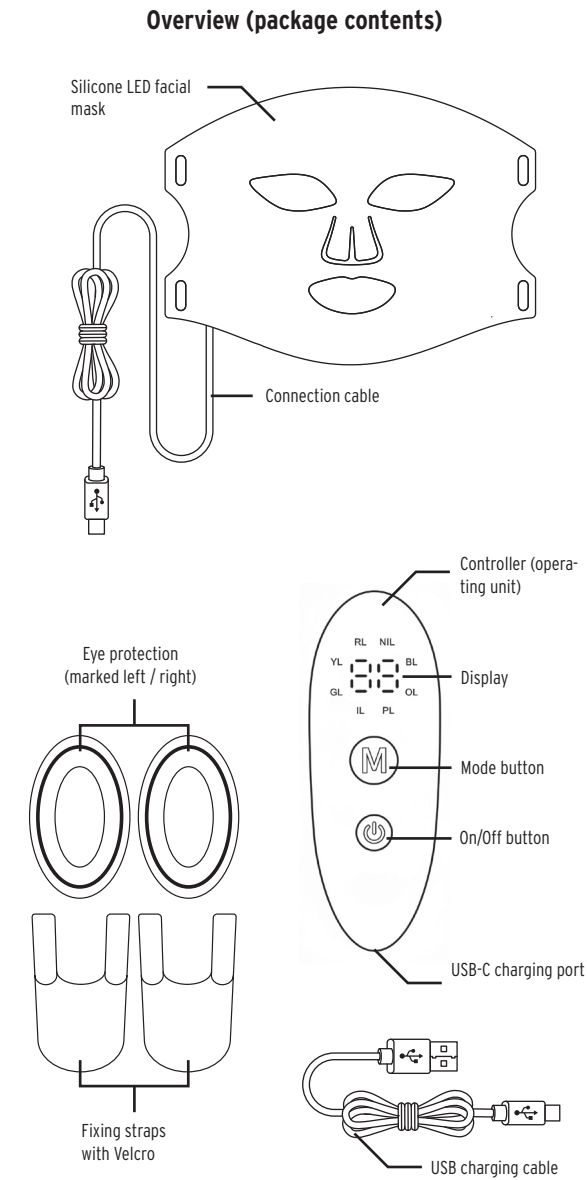EG - DECLARATION OF CONFORMITY

The product described in this instruction manual complies with all relevant harmonized requirements of the EU. The complete declaration of conformity and instruction manual can be found online at www.jaxmotech.de/downloads.

About these instructions

This product is fitted with safety features. Please read the safety instructions carefully and use the product only as described in this manual to avoid accidental injury or damage. Please keep these instructions for future reference. If you pass this product on to someone else, please include these instructions as well. Symbols used in these instructions:

 This symbol warns you of a risk of injury.

 This symbol warns you of the risk of injury from electricity.

Key terms in these instructions:
DANGER: Indicates an imminent risk of serious injury or death.
CAUTION: Risk of minor injury.
NOTE: Warns of possible damage to property.

 Additional information is marked with this symbol.

Safety instructions

Purpose

This LED facial mask is designed for use on healthy, uninjured facial skin. It uses various LED light colours to support skin care (e.g. to improve skin appearance or well-being).

Note: This product is not a **medical device** as defined by the EU Medical Device Regulation (MDR) and does not replace **medical treatment**.

Not suitable for:

- Children and young people under 18 years
- Individuals with known light sensitivity
- Individuals with active skin conditions, open wounds, or severe inflammation in the facial area
- Pregnant or breastfeeding individuals (use only after consulting a doctor)

Any other use, misuse or improper use may result in damage to the appliance or personal injury.

Danger to children and people with reduced ability to operate equipment.

- This appliance may be used by persons aged 18 and over with reduced physical, sensory or mental capabilities, or a lack of experience and/or knowledge, provided they are supervised or have been instructed in the safe use of the product and understand the resulting risks.
- Cleaning and maintenance by the user must not be carried out by children.
- Children must not play with the appliance.
- Keep children away from the appliance and the mains lead.
- Always keep the appliance and all accessories out of the reach of children.
- Before carrying out any cleaning, unplug the appliance, allow it to cool down and follow the instructions in the ‘Cleaning and Care’ section. Do not immerse the product in water when cleaning it.
- People who are not using the appliance (in particular children) must not be present whilst it is in use.
- The device is not intended for outdoor use.
- Only use the appliance in dry indoor areas.
- Do not use the appliance in environments where there is a risk of explosion or fire.
- Never look directly into the LED light.
- High brightness. Do not look into the spotlight. This can result in eye injuries.
- This product is intended for domestic use only.
- Eye or skin irritation, or even injury, may occur following prolonged exposure.
- Never leave the product unattended whilst it is switched on.
- Never immerse the electrical parts of the product in water or other liquids during cleaning or operation. Never hold the product under running water.

- Maintenance and repairs must be carried out by an authorised service centre.

Information about rechargeable batteries

- This appliance contains built-in rechargeable batteries.
- Before disposing of the appliance, contact an electrician or your retailer to have the batteries removed and disposed of in an environmentally responsible manner.
- Disconnect the appliance from the mains before removing the batteries.
- Never expose batteries to excessive heat (e.g. direct sunlight, fire). The batteries could explode.
- Never short-circuit the battery terminals.
- Do not throw the battery or the appliance into a fire. Risk of explosion and fire!
- Never leave the appliance's battery charging unattended.

Usage

 **DANGER** to children

Keep children away from packaging materials. Risk of suffocation!

Before first use

- Remove all packaging material.

Light modes and effects

Mode	Colour	Effect
RL	Red	Promotes collagen regeneration, improves elasticity, reduces wrinkles and fine lines
NIL	Near-infrared	Promotes circulation
BL	Blue	Refines pores, supports skin balance
OL	Orange	Soothes the skin, assists with pigmentation issues
PL	Violet	Lightens dark spots, promotes collagen
IL	Indigo	Supports skin metabolism
GL	Green	Reduces blackheads and pigmentation spots
YL	Yellow	Improves rough skin, gentle collagen promotion

Charging the controller

- Do not use the device while it is charging.
- Flashing indicator = charging in progress
- All indicators illuminate continuously = fully charged
- Fully charge before first use.

Assembling the mask

Attaching the eye protection

1. Place the mask with the front side facing down.
2. Insert the eye protection with L on the left eye and R on the right eye.
3. The letters must face outward.

Attaching the straps

1. The straps are universal (no left/right).
2. The soft flannel side faces inwards.
3. Thread the straps through the side openings and secure with Velcro.
4. Once the flannel straps are secured, the large Velcro strips can be closed around the head.

Connection

Connect the face mask to the controller.

OPERATING THE CONTROLLER

Button functions

- **On/Off & time selection:**
 - Press and hold: On/Off switch
 - Press briefly: Time selection (**5 minutes or 10 minutes**)
- **Mode button**
 - Switch between the 8 light modes by pressing briefly (RL → NIL → BL → OL → PL → IL → GL → YL)

Display:

- Illuminated mode indicator = mode enabled
- Flashing indicator = low battery
- Automatic switch-off

First use – skin compatibility test

It is recommended to conduct a test before first use:

1. Select a light mode.
2. Test the device on a small area of skin (e.g. forearm) for 5-10 minutes.
3. If no skin reactions occur within 24 hours, the product can be used on the face.
4. If in doubt, consult a dermatologist before further use.

Application steps

1. Cleanse the face thoroughly and dry it.
2. Put on the mask and adjust the straps.
3. Press and hold the controller to switch the device on.
4. Select the light mode and set the application duration (**5 or 10 minutes**).
5. The device will automatically switch off after the time has elapsed.
6. Remove the mask and apply your usual skincare.

Recommended usage:

- 3-5 applications per week
- Use for a maximum of 10 minutes per session

Do not use in cases of:

- Pregnancy or breastfeeding
- Light sensitivity or taking photosensitising medications
- Eye or skin conditions (e.g. lupus, photosensitive eczema, albinism)
- If you experience burning or discomfort, discontinue use immediately.

Cleaning and storage

1. Clean the appliance regularly with a dry, soft cloth.
2. Do not use aggressive cleaning products.
3. Keep accessories clean and dry after use.

Storage

- Store the appliance out of the reach of children.

Technical specifications	
Model	: FM 008
Item number	: 763 579
USB input	: 5 V === 1 A
Number of LEDs	: 280 LEDs
Maximum power	: 4 W
Automatic runtime	: 5 and 10 minutes
Battery: Lithium-ion	: 3.7 V === 1500 mAh
Watt-hour rating	: 5.55 Wh
Charging time	: Approx. 3.5 hours
Operating time	: max. 80 minutes
Protection class:	III
Information about the built-in battery:	
Manufacturer	: Shenzhen Julang Energy Industry Co., Ltd 309, Building A3, No. 19, Youmagang Industrial Zone, Shiwei Community, Matian Street, Guangming District, Shenzhen City, China
Email	: julanglian0213@163.com
Battery model	: JL 18650 1500mAh
Date of manufacture	: 2026/06 Made in China
Manufactured by	: Jaxmotech GmbH Ostring 60, 66740 Saarlouis, Germany info@jaxmotech.de

As part of product improvements, we reserve the right to make technical and visual changes to the item.

Disposal

The product, its packaging and the built-in rechargeable battery have been manufactured from valuable materials that can be recycled. Recycling reduces the amount of refuse and preserves the environment. Dispose of the **packaging** at a recycling point that sorts materials by type. Make use of the local facilities provided for collecting paper, cardboard and lightweight packaging.

Appliances, single-use batteries and rechargeable batteries marked with this symbol must not be disposed of along with household waste! You are legally obliged to dispose of old appliances separately from household waste. Electronic appliances contain hazardous substances and, if stored or disposed of improperly, may cause harm to health and the environment. Information about collection points where old appliances can be disposed of free of charge is available from your local authorities.

Flat single-use and rechargeable batteries must be handed in at a local authority collection point or returned to a battery retailer. Tape off the contacts of lithium single-use/rechargeable batteries prior to disposal.


Warning! This appliance contains rechargeable batteries that, for safety reasons, are firmly built into the appliance and cannot be removed without destroying the casing. Improper disassembly leads to safety hazards. For this reason, take the appliance intact to a collection point where the appliance and battery will be disposed of properly.

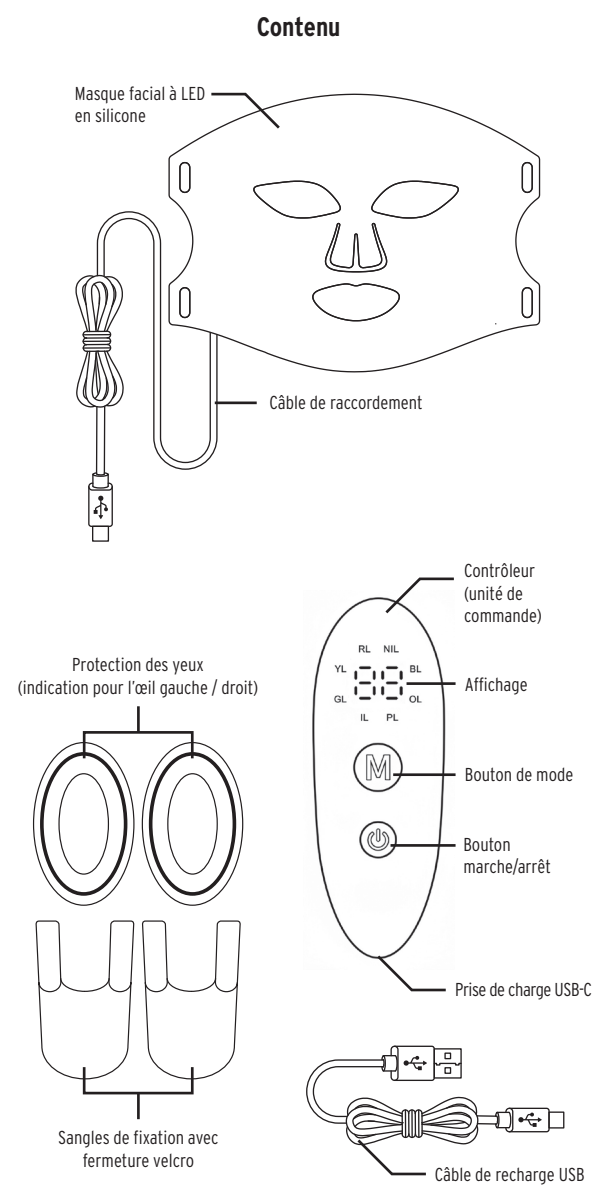Model	: FM 008
Manufactured by	: Jaxmotech GmbH, Ostring 60, 66740 Saarlouis, Germany info@jaxmotech.de

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Hersteller : Shenzhen Julang Energy Industry Co., Ltd
Hergestellt von : Jaxmotech GmbH
Hergestellt von: Jaxmotech GmbH,